HOW TO USE BACK AND NECK HEATING PAD

DESCRIPTION

- 1) OFF button
- 2) Position 1 – 1st level of temperature
- 3) Position 2– 2nd level of temperature
- 4) Position 3– 3rd level of temperature

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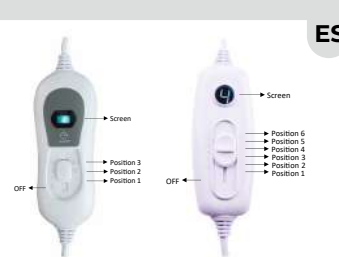
FR

OPERATING

- Put on the electric heated vest using the tether strap and Velcro for the neck.
- Connect the plug to the socket at the bottom of the heating pad and the other end to an electricity supply (AC220-240V).
- To pre-heat the heating pad, select the third setting for 15 minutes. The temperature will gradually rise.
- After pre-heating, turn the control button to the heat setting you require.

Note: The Heating Pad is set to run continuously for 90 minutes. If you wish to extend that time, it is advisable to let it cool down for 5 minutes before using it again.

USER MANUAL

1. Do not use by children under the age of three (from 0 to 3 years).
2. Do not use the electric heating pad folded or wrinkled.
3. Do not use if wet.
4. Electric and magnetic fields from this electric heating pad may interfere with the operation of pacemakers.
5. You should therefore not cover switch or place it on the heating pad when the latter is used.
6. Do not pull, twist or bend the cables.
7. This heating pad must be regularly checked for the presence of any sign of wear or damage.

CLEANING AND MAINTENANCE

When you are not using the electric heating pad and before cleaning it, follow the instructions below:

- Make sure the control button is in the "OFF" position.
- Disconnect the power supply.
- Disconnect the plug linking the control to the device.
- Let the heating pad cool down before cleaning.

Small stains can be cleaned with a cloth or a damp sponge and, if necessary, with liquid detergent for delicate fabrics.

Warning: the electric heating pad should not be dry-cleaned, bent, dried in a dryer or machinepressed or ironed.

If the heating pad is very dirty, it can be washed in a washing machine. Set it to the 30°C program for delicate fabrics.

Use a detergent for delicate fabrics and measure it out according to the manufacturer's instructions.

The heating pad should not be washed in the washing machine more than 5 times during its life.

Never switch the electric heated vest on to dry it.

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APPLIANCE CO.,LTD
Imported by/Importado por: Grupo R. Queraltó S.A.