Manual de instrucciones

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Massage Gun

PS-01

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User manual in English **EN**

ÍNDICE

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INDEX

1. Parts included	8
2. Product features	9
3. Operating instructions	10
4. Battery instructions	11
5. Warning and precautions	11
6. Maintenance	12


This item is not a medical device and is not intended for the treatment or prevention of disease.

1. Included parts


1 TYPE A HEAD

Suitable for different areas of the body.

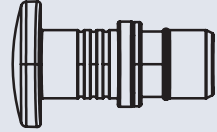
2 TYPE B HEAD

Suitable for hands, among other areas of the body.

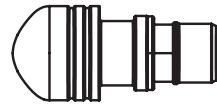
3 TYPE C HEAD

Suitable for different areas of the body.

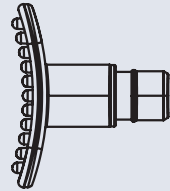
4 TYPE D HEAD

Suitable for different areas of the body.

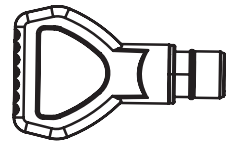
5 TYPE E HEAD

Suitable for arms, waist, thighs and cervical area.

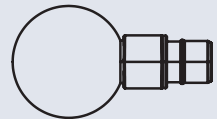
6 TYPE Y HEAD

Suitable for the cervical area and ankles.

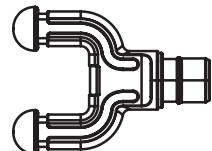2. Product features

Name	
1	Charing socket
2	Power Switch
3	Speed up
4	Speed Down
5	Batery indicator
6	Speed level

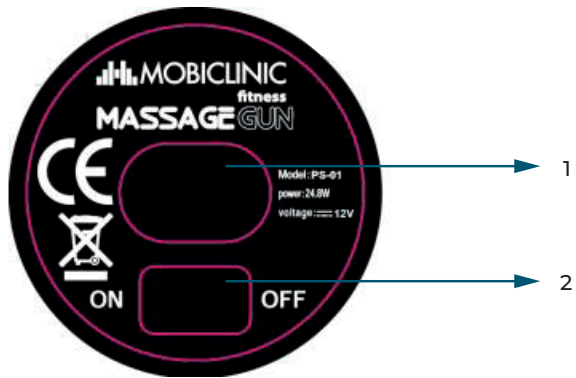3. Operating Instructions

1 Charging socket

Please charge the battery fully when using this device for the first time. Use the original charger.

Warning: Make sure that the device is not running during charging.

2 Switch

To turn the device on, turn the switch to the left (ON), and to turn the device off, turn the switch to the right (OFF).

3 Speed Up

Press the + symbol on the display to increase the massage speed level.

4 Speed Down

Press the - symbol on the display to decrease the massage speed level.

5 Battery indicator

Displays the charge level of the device.

6 Speed Level Indicator

Starting from 1350rpm/min at level 1, each level increases by 50rpm/min up to level 30 with a maximum of 2800rpm/min.

4. Battery Instructions

1. This device is equipped with 12.5V/1500mAH lithium battery, which is safe and reliable, and can work continuously for 2.5 hours at the lowest speed level position.
2. To charge the battery fully it should be on charge for 3 hours uninterruptedly.

Warning: The battery of this appliance is not removable. If you need to replace the battery, contact a professional technician.

5. Warning and precautions

Before using this device, be sure to read these warnings and cautions carefully:

- Do not use this product for infants, minors, pregnant women, persons wearing pacemakers or having metallic contents in the body.
- Keep the product away from liquids, heat sources and sources of ignition.
- It is prohibited to perform any disassembly and modification on the product.
- The device cannot be used on the neck and head.
- The time of continuous massage of the same part must not exceed 60 seconds.
- Be sure to install a head before using the device. According to the different parts of the body, choose the appropriate massage head.
- Keep fingers, hair or any part of the body away from the shaft and the back of the massage head.
- Do not immerse the product in water or allow water to enter the product's vents.
- Do not discard or misuse the device.
- Never leave the product unattended.

6. Maintenance

1. Clean the surface of the unit with a slightly damp cloth and wipe it dry with a dry cloth.

Made in P.R.C

ZHETJIANG RONGSHUN
TECHNOLOGY CO., LTD
No. 399 Zhengnan North St,
High-Tech Industrial Park,
Jinyi New Urban District,
Jinhua, Zhejiang, China

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GRUPO R. QUERALTÓ S.A.